

**Sitting
& Standing
Support!**

Stability Where It Matters Most

Wide Base of Hand Support: Enhances balance and postural control during transitions by providing stability at every phase of movement.

Support Directly In Front: Pivoting bar delivers assistance exactly where clients need it – ideal for sit-to-stand, pivot, and bedside transfers.

Solid Support for ADLs: Stable enough for leaning against while dressing, managing clothing, personal hygiene, and toileting tasks.

Reduces Caregiver Strain: Enables clients to reposition and stand with minimal hands-on help. Decreases physical demands on staff and lowers risk of caregiver injury.

Tool-Free, Floor-to-Ceiling Installation: No drilling required – pressure-fit design installs quickly, even in rental homes or short-term rehab environments.

Engineered for Stability: Supports up to 300lbs and fits ceiling heights from 7–10 feet. Install anywhere—beside beds, toilets, chairs, or therapy stations.



Every Surface a Secure Grip

Clean, obstruction-free design – no jack screws or mechanical protrusions – ensuring a confident grip anywhere.



8 Locking Support Positions

The bar locks every 45°, providing consistent and safe support from multiple angles. Built with the same SuperBar technology, trusted for over 30 years.



Effortless Pivoting Bar


Easy to pivot with just one hand. Designed for users with limited hand strength or dexterity – no buttons, levers or squeezing required.

uPole™ + Pivot Bar

1 Turn
jackscrew



2 Press top plate
to ceiling

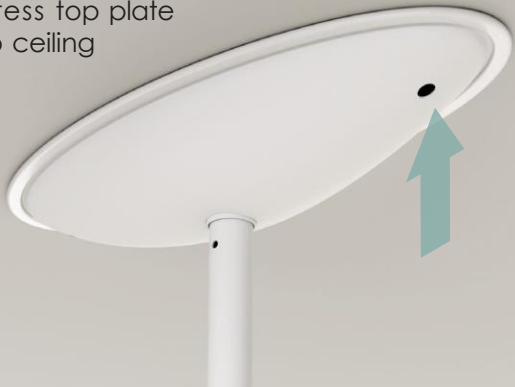

3 Lock with jamnut



4 Conceal with cover

